

# K-scope Head Light with Camera

(Model : KS-07)



All-in-One

Cordless

Rechargeable



# K-scope

## Head Light with Camera

(Model : KS-07)

01

### It has “all-in-one” style

Since the camera is built-in, there no need to equip/install another separately

02

### It has Cordless

Able to record, without power, in MicroSD card and rechargeable battery cordlessly.

03

### It has dual usage (battery and/or power)

It has rechargeable battery. You may use cordlessly. When in emergency or longer usage, it is also capable to connect power source. (option)

04

### It has LED lamp

It is built-in LED lamp(5500K), it lights brighter and clearer light for over 50,000 hours.

05

### 1.7X Loupe is provided as free accessory

Simple installable free 1.7X Loupe is provide.

06

### It is able to record and/or demonstration

For several purposes and usages, it is able to record using MicroSD card. Also able to demonstrate using video cable.

07

### K-scope Viewer Program is provide

It is provided own-use “K-scope Viewer” program. This is capable to make “time-reset” and “voice cancellation”



Demonstrate using video cable



## Convenient Functions


1. Sound (and/or alarm)
  - ① When in no MicroSD card
  - ② When in low battery voltage
  - ③ When in camera on/off
2. Adjust camera focus
3. Adjust camera angle
4. Automatic time-resetting (K-scope Viewer)
5. Voice cancellation (K-scope Viewer)
6. Extension video Cable (Option)

## Usage

1. For research and/or presentation of theses
2. For field-training by demonstration
3. When in dispute (blackbox concept)
4. For insurance proof



Light and Camera

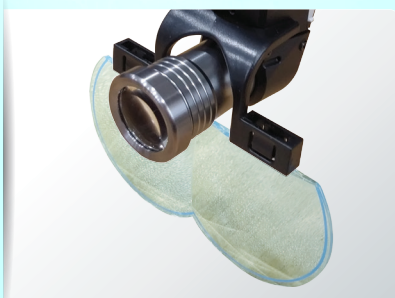

1.7X Loupe



"K-scope Viewer" Program



## Components



① Headlight main body	1
② 1.7X Loupe	1
③ Loupe Holder	1
④ Rechargeable Battery	2
⑤ Battery charging holder	1
⑥ Charging Adapter	1
⑦ MicroSD Card (included "K-scope Viewer")	1
⑧ MicroSD Card reader	1
⑨ Video Cable	1
⑩ USB Cable	1
⑪ Manual	1

## Optional Accessories



- ① Video and power Cable
- ② Extension Video Cable (3M)

## Specification

- Lamp : 5500K Power LED
- Battery : 2,600mA Li-ion Battery
- Camera
  - Sensor : 1/4" CMOS WXGA (1 megapixel) HD sensor
  - Memory : MicroSD (SDHC 2GB~64GB)
  - Resolution and Frame : 1280\*720 HD / 30fps
  - Output : TV OUT (composite)
  - Sound : Built-in Mic.
- Lens : Scientifically arranged three optical lens

